



ADD-ON DRIVE  
COMPATIBLE



### Where stability meets style!

Invacare **Action<sup>4</sup><sub>NG</sub>** offers a flexible range of models to accommodate different customer needs and is particularly suitable for taller or heavier people.

With a frame designed for improved stability and comfort, the stylish **Action<sup>4</sup><sub>NG</sub>** offers the same cross compatibility of parts with other models in the **Action<sup>4</sup><sub>NG</sub>** family allowing customers to easily maintain and update their wheelchairs when required.



**Yes, you can.®**



### Easy rolling and maneuverability

The frame design of **Action<sup>4</sup><sub>NG</sub>** helps improve stability resulting in easier wheeling and manoeuvring by maximising the customer's position in the chair.

It is easy to adapt the **Action<sup>4</sup><sub>NG</sub>** to other types of wheeling characteristics to suit various needs by changing the rear wheel position and altering the seat position. It features optional seat adjustments to aid postural support for day long use.

### Compatibility with other Invacare wheelchairs

Almost all parts of the **Action<sup>4</sup><sub>NG</sub>** frame are interchangeable with other models of the **Action** family and even with other Invacare wheelchairs, including *Rea*™. All parts such as armrests, legrests, wheels etc, are compatible with each other. The **Action<sup>4</sup><sub>NG</sub>** is also compatible with most of the **Alber** product range - please see the **Alber** brochures for further details.

### Easy transport with low weight

By removing the rear wheels and folding the frame, the **Action<sup>4</sup><sub>NG</sub>** easily slides in behind the front seat or slips into the boot of a car. Without rear wheels, the **Action<sup>4</sup><sub>NG</sub>** weighs 8.3 kg. Simple handling and folding makes for great ease of transport.

Even the most basic **Action<sup>4</sup><sub>NG</sub>** model covers customer needs, as almost everything, from the height of

the backrest, armrest and footrest is adjustable.

From the **Action<sup>4</sup><sub>NG</sub>** Tall (125 kg) Transit model for use with an attendant, to the self-propelling model or the groundbreaking **Action<sup>4</sup><sub>NG</sub> Heavy Duty** version (with a user weight of 160 kg) the **Action<sup>4</sup><sub>NG</sub>** offers total flexibility in a chair.

## Features and Options



Available as either a standard Transit version or Transit HD.



### Adjustable backrest upholstery

As an option, the backrest can be equipped with tension-adjustable straps, shaping the backrest to individual needs.



### Seat-depth adjustability

The seat depth has up to 5cm adjustment with the upholstery.



### Backrest bar

With tension adjustable upholstery a backrest bar is available for extra rigidity.

## Features and Options



### Hangers and legrests

Hangers come as angle-adjustable, fixed 80 degree angle and amputee version. Use depth-adjustable, fixed or one piece footrest with your hanger.



### Optional rear wheel positions

New active position plate offers better balancing possibilities for superior driving performance.



### Headrest

Can be adjusted on the backbrace in width and depth to provide optimum positioning.



### Hub brakes

Hub brakes available as an option, helping assistants to stop the wheelchair and reduce speed.



### Anti-tip device

Prevents wheelchair from tipping over and is height adjustable to suit the different seat heights available.



### Amputee support

The amputee support is both height and angle adjustable.



### Heavy Duty Version 160 kg

A wider choice of seat widths (50.5, 55.5 & 60.5cm), reinforced chassis and a greater carrying capacity (maximum user weight 160 kg).



### Hemi version

Dual handrim propulsion and hemi armrest are available.


## Technical data

Invacare <b>Action<sup>4</sup>NG</b>	38 / 40.5 / 43 / 45.5 / 48 cm	40 / 45 / 50 cm	41 / 43.5 / 46 / 48.5 / 51 cm	38 - 46 / 39 - 51 cm*	20 - 23 - 26 / 19 - 29.5 cm	35 - 48 cm
Invacare <b>Action<sup>4</sup>NG HD</b>	50.5 / 55.5 / 60.5 cm	40 / 45 / 50 cm	41 / 43.5 / 46 / 48.5 / 51 cm	38 - 46 / 39 - 51 cm*	20 - 23 - 26 - 28 cm / 19 - 29.5 cm	35 - 48 cm
Invacare <b>Action<sup>4</sup>NG</b>	0° - 80°	7° / 0° - 30°**	Seat width + 20 cm	90 - 103 cm	105 - 111 cm	83 - 90 cm
Invacare <b>Action<sup>4</sup>NG HD</b>	0° - 80°	7° / 0° - 30°**	Seat width + 20 cm	90 - 103 cm	105 - 111 cm	83 - 90 cm
Invacare <b>Action<sup>4</sup>NG</b>	30.5 cm	15.2 kg	125 kg	8.3 kg	Crash tested ISO 7176-19**	
Invacare <b>Action<sup>4</sup>NG HD</b>	30.5 cm	17 kg	160 kg**	8.3 kg	Crash tested ISO 7176-19**	




\* Reclining backrest  
 \*\* Our wheelchairs comply above ISO standards and with ISO 10542, and have additional features providing increased levels of occupant security and safety whilst travelling in a motor vehicle. Wheelchair users should however transfer to the vehicle seat and use the vehicle installed restraint system whenever it is feasible.

## Frame colours

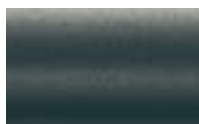

Slate Grey

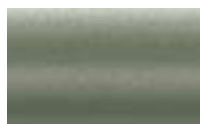

Ice Brown



Ocean Blue



Anthracite Black



Dark Red

## Upholstery



Black Nylon

INVACARE Limited  
 Pencoed Technology Park - Pencoed -  
 Bridgend CF35 5AQ  
 Tel: +44 1 656 776 222  
 Fax: +44 1656 776 220  
 E-mail: uk@invacare.com  
 www.invacare.co.uk  
 Sales Order E-mail:  
 ordersuk@invacare.com

INVACARE Ireland Ltd  
 Unit 5 - Seatown Business Campus  
 Seatown Road - SWORDS - County  
 Dublin  
 Tel. +353 1 8107084  
 Fax +353 1 8107085  
 Email: Ireland@invacare.com  
 www.invacare.ie

© 2011 Invacare International Sàrl  
 All rights reserved.

All information quoted is believed to be correct at time of print.  
 Invacare reserves the right to alter product specifications without prior consultation.

Action4NG - UK - 07/2011  
 1528015



Yes, you can.®



Please notice colors may vary slightly from those displayed above.