

# Performance, style and comfort

**S-Ergo 125**



The Ergo 125 is a lightweight aircraft grade aluminium wheelchair which breaks down to 9kg (20lbs). The lightweight frame is braced by a system of double cross braces which increase the wheelchair's strength and durability. The streamline oval tubes give the wheelchair its unique look. The Ergo 125 is packed full of innovative features that you would expect from a Karma wheelchair including, flip back armrests, slowing brakes, quick release wheels, detachable upholstery and the patented Ergo seating.



▲ Self propel



▲ Transit



Pressure redistribution



Reduced sliding



**TALL**  
**S-Ergo 125**



Extended seat depth to assist taller users in need of mobility.

## HOW DOES IT WORK?




**TRADITIONAL SEAT**

High risk of decubitus ulcers could be occurred because pressure points are concentrated on hipbone area.

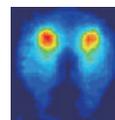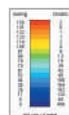

**S-Ergo SEAT**

The S Ergo system promotes pressure relief by spreading weight over a greater area. This also increases stabilization and reduces sliding.

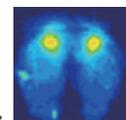



The streamline oval tubes not only make the wheelchair look great but also add strength to the wheelchair structure.

## FEATURES

### Swing away footrest



Swing in/out and detachable footrest adds convenience when the user transfers in tight areas.

### Flip-back armrest



Flip-back armrest offers greater convenience if the user side transfers.

### Quick release rear wheels



### Outrigged armrests



The outrigged armrests are designed to provide a comfortable arm support as well as maximize seating area.

### Ergo II Armrest pad



The smooth edge of arm pad provides comfortable support but does not restrict the user from propelling the wheelchair.

### Ergo handrim

Patents awarded in multiple territories.



The innovative Ergo handrim offers handling comfort and assists with wheelchair propulsion.



S-Ergo entire series have seat upholstery and pads treated with AEGIS microbe shield to reduce stains and odour.



Foldable and extremely compact

## OPTIONS



Elevating legrest



Stump support



Foldable headrest



Spoke guard



Travel bag

Model	items Width	F/R wheels	Options										Max. User Weight
			Legrest	Stump	Headrest	Spoke	Travel	Seat	Back	Wheels	Armrest	Handrim	
S-Ergo 125	16"	7"-24"	40	43	57	30	21	91	50	41	91	12.9	115kg
		7"-14"	40	43	59	31	21	99	50	41	91	14.3	115kg
	18"	7"-24"	45	43	62	30	21	91	50	41	91	12.4	115kg
		7"-14"	45	43	63	31	21	99	50	41	91	13.1	115kg
	20"	7"-24"	50	43	67	30	21	91	50	41	91	14.5	115kg
		7"-14"	50	43	68	31	21	99	50	41	91	14.7	115kg
S-Ergo 125 Tall	18"	7"-24"	45.7	48	64	32.5	21	107	50	41	90.5	14.6	115kg
		7"-14"	45.7	48	62	32.5	21	96	50	41	90.5	13.5	115kg
	20"	7"-24"	50.7	48	69	32.5	21	107	50	41	90.5	14.6	115kg
		7"-14"	50.7	48	67	32.5	21	96	50	41	90.5	13.5	115kg

±1cm ±0.5kg

KARMA reserves the right to modify information herein without further notice. - 104-405A2 Rev.03 Jun.2018