

We continue to expand our rollator and entire product portfolio

For several years now, TOPRO has continually updated, or added to, its rollator line-up – e. g. the Olympos ATR (All-Terrain-Rollator) or Pegasus, our lightweight carbon model. This time around, we are proud to present you with the TOPRO Hestia, our all-new indoor-rollator with never-before-seen features and benefits (for more details, turn to pages 6 to 10).

As the market leader for pioneering mobility solutions we don't stop there: our product portfolio now includes electric mobility scooters as well as AssiStep, a truly innovative, mechanical "rollator" for stairs (stair walker) – please refer to our website or the respective brochures for more information.

Developing and producing quality mobility aids has been our focus for over 30 years. We take pride in what we do and are happy to see how our products improve our customers' daily lives in the here and now.


Looking to the future, beyond mere product benefits, the environment also plays a big role in our thought process. with sustainability playing a central role: Our factory in Norway uses hydropower and environmentally friendly aluminium to minimize the carbon footprint. Our injection-moulding plant uses clean materials that are recyclable to contribute to the circular economy. We have also implemented a "zero

waste"-principle to ensure no plastic contaminates the natural environment for the benefit of future generations.

Our goal, through our products, is to make your life easier. Please stay healthy, active – and "keep on moving"!



Asbjørn Eskild **CEO TOPRO Group**

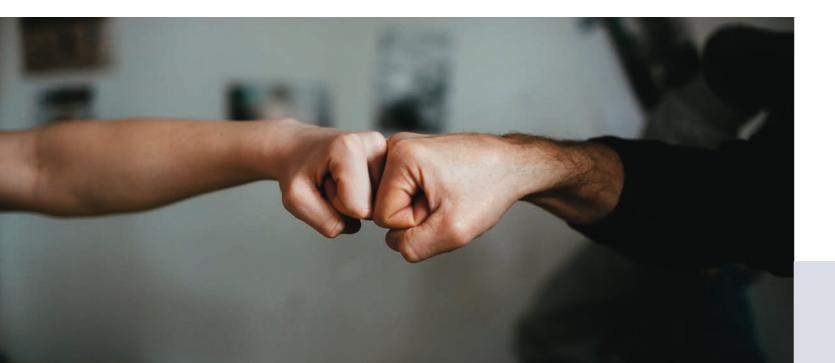

Advantages of the TOPRO rollator family

Whether it's product features or design – we have continually improved our rollators for over 30 years.

We are convinced our products provide you with an added value – by combining smart functionality and modern design with high quality and safety.

Unique memory function

Saves the preferred preset handle height, for easy re-use. (Standard on Troja 5G, Olympos ATR and Troja Original)



Ergonomic handles

The TOPRO Ergo Grip ensures correct walking posture and comfortable armrests while seated. Together with the soft tyres it provides double shock absorbing effect for stable driving comfort. (Standard on TOPRO Pegasus, Troja 5G, Olympos ATR and Odyssé)

Easy to fold

Simply pull the strap on the seat with only one hand. Practical to transport, stays firmly on its wheels and takes up little space when stowed.

Easy wheel change

The interior braking system (TOPRO IBS) in combination with the TOPRO Quick Release system for easy wheel change, allows the use of different tyre treads. (Standard on TOPRO Pegasus, Troja 5G, Olympos ATR and Odyssé)

Tilt function and edge guard

Enables the user to easily negotiate thresholds and kerbs, while keeping both hands on the rollator for safety.

Find the rollator that fits you best:

Our rollator guide assists you with selecting the right rollator for your individual needs:



TOPRO PEGASUS





Modern, slender and lightweight carbon rollator with functionality and elegant design. This design combines to guarantee the highest level of quality, safety and stability at a weight of just 6.1 kg*.

Tried-and-tested features have been further developed:

- User-friendly method for setting the height of the handles using a numeric scale.
- Partly integrated brake cables and modern reflective panels on the frame.
- Elegant detachable shopping bag, which, with its retracted position provides increased forward view to improve user safety.

- Easy handling and manoeuvrability.
- Reflectors on the frame, lock wheels and shopping bag provide 360 degree visibility for safety in the dark.
- Shock-absorbing tyres and ball bearings in wheels and forks for a comfortable ride.

Measurements, colour options and other product details are available on page 18.

Available in colour:





Steel blue



Wine red - NEW!









* Weight including shopping bag 6.6 kg.









Stylish, easy-to-fold space-saving indoor rollator – offering you that extra bit of self-reliance and comfort at home. With a slew of clever features, the Hestia enables you to move about safely, get up and sit down more easily as well as to tranpsort meals and other items securely.

You can look forward to these innovative features:

- Integrated, extendable rise-up handles provide reliable support for getting up and sitting down. The handles retract for unobstructed walking.
- Multi-purpose tray can easily be lifted and carried with one hand with its retractable handle, to safely transfer food and drinks.
- The slidable tray can be moved to the walker side, making it easier to reach
- A cup holder and a small groove for your smartphone or tablet come standard. Clever: If you place your loudspeakers in the holder, their volume gets amplified!
- Handrail with a continuous brake lever allows for variable grip positions and safe handling while manoeuvering the compact Hestia.





Measurements, colour options and other product details are available on page 19.













Hestia: The new benchmark for indoor rollators


The TOPRO Hestia is ready to go – no assembly required It means you can get moving and use these practical details right away:



TOPRO Ergo Handle – Ergonomic and height adjustable, the one-piece handle-bar allows variable grip positions as well as single-handed operation.



Continuos Brake Lever – Safe and simple braking from different grip positions. Integrated parking brake.

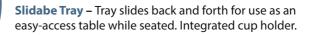

Carry with one Hand – Tray is easily detachable and balances itself when carried with one hand.



Basket – Discreet storage for personal items underneath tray.





Tilt-function – To overcome small obstacles (e. g. thresholds)



TOPRO Hestia is both a sturdy and compact, lightweight indoor rollator, easy-to-manoeuvre around your home. When you don't need it, it folds for compact storage (including the tray) and stands on its own, with a small footprint.

TROJA 5G



Handling comfort and stability – in true TOPRO-style. The TOPRO Troja 5G sets itself apart with its modern, slender design that goes hand-in-hand with a high level of quality and a secure driving experience.

It's the little things that make the difference – for improved usability and comfort:

- User friendly method for setting the height of the handles using a **numeric**
- Shock-absorbing tyres provides a comfortable ride and smooth suspension.
- Bell to attract attention when necessary. Easy and safe to operate from the handle.
- Easy handling and manoeuvrability due to ball bearings in forks and wheels.

- Detachable shopping bag with reflective band, zipper and shoulder strap. Closeable and holds up to 10 kg. Two inner pockets for small items.
- Easy to fold, stays firmly on its wheels and takes up little space when stored.

Measurements, colour options and other product details are available on page 20.

Model M & S available in on-trend matt coated colours:













OLYMPOS ATR

THE ALL-TERRAIN ROLLATOR



With large front wheels the Olympos ATR provides reliable and comfortable support on uneven roads and cobblestones as well as on sand or gravel foundations. The perfect choice for outdoor activities.

Innovative features to provide maximum support for outdoor use:

- Wider design to ensure stability in an outdoor environment.
- Large and shock-absorbing tyres to enable a comfortable ride and maintain stable and smooth manoeuvrability on rough terrain.
- Detachable spacious water-resistant **shopping bag** in on-trend grey fabric print. With reflective band and zipper. Closeable and holds up to 15 kg.
- User friendly method for setting the height of the handles with numeric scale.
- Easy to ride and manoeuvre due to ball bearings in forks and wheels.

Measurements, colour options and other product details are available on page 22.

Model M & S available in on-trend matt coated colours:















TROJA ORIGINAL





TOPRO's basic rollator with state-of-the-art design and high level of quality and features. The slender and modern design is infused with the functionality needed for a safe walk and an active everyday life.

The functionality is based on TOPROs high and continuous quality over many years:

- Lightweight at only 6.8 kg, and with a weight bearing capacity of 165 kg.
- Innovative sideways folding techno**logy** that is easy to operate with one hand by pulling the strap on the seat.
- Ball bearings in forks and wheels to ensure an easy ride and manoeuvrability.
- Practical net basket for easy access to items. Reflective band for safety and visibility in the dark.

- Available in three different sizes.
- Valued by users for being easy to handle and operate, having great manoeuvrability which makes it a trusty companion in everyday situations.
- Numeric scale for easy height adjustment of the handles.

Model M, S and X available in:





Wine red

Measurements, colour options and

other product details are available on page 21.









Accessories



Crutch holder can be fitted on either side.

Troja Original / 5G, Olympos ATR, Hestia – grey Art.No. 815358 TOPRO Pegasus Art.No. 815236 Olympos Art.No. 814040 Odyssé Art.No. 814457

TOPRO Pegasus, Troja 5G, Olympos ATR, Troja Original, Odyssé, Troja Neuro, Hestia



Snaps in securely, easy to remove and clean, carries up to 5 kg

TOPRO Pegasus, Troja Original / 5G, Odyssé Art.No. 814728 Olympos ATR M Art.No. 815212

TOPRO Pegasus, Troja 5G, Olympos ATR, Troja Original, Odyssé



Back support with padding ergonomically height adjustable according to preferred seating position. Available in two lengths.

Short 67 cm Art.No. 815270 Art.No. 815272 Long 75 cm

TOPRO Pegasus, Troja 5G, Olympos ATR, Odyssé, Troja Neuro



LED lamp with clip on and USB **plug** for safety in the dark. Chargeable.

Art.No. 815260

TOPRO Pegasus, Troja 5G, Olympos ATR, Troja Original, Odyssé, Troja Neuro



Back support with padding ergonomically height adjustable according to preferred seating position. Available in two lengths.

Art.No. 815271 Short 68 cm Art.No. 815273 Long 76 cm

Troja Original



Umbrella with attachment keep dry under the umbrella in rainy weather or use it as a parasol.

Black Art.No. 814068 Red Art.No. 814099

TOPRO Pegasus, Troja 5G, Olympos ATR, Troja Original, Odyssé, Troja Neuro

Additional accessories can be found on our website.





Technical data for TOPRO Pegasus



| NEW | VEAR |
|-----|---------|
| | 4ARRANT |

Available in colour:

Wine red - NEW!

Technical data for Hestia



| ight 15 | 50–200 cm |
|---|----------------------------|
| user | 130 kg |
| bag | 10 kg / 15 l |
| ight 7 | 78–97.5 cm |
| dles Ergonomic, serveral holding (TOPRO | g positions) ErgoGrip) |
| rake | 7 cm |
| dles | 47 cm |
| idth | 43.5 cm |
| ight | 63 cm |
| ircle | 84 cm |
| ight | 101 cm |
| idth | 64 cm |
| ngth | 65 cm |
| eels | 20 cm |
| eels | 3.6 cm |
| lded | 65 cm |
| lded | 27 cm |
| lded | 80.5 cm |
| ×h) 65 × 27.3 | 3 × 80.3 cm |
| bag) | 6.1 kg |
| bag) | 6.6 kg |
| fuse Indoors 8 | & outdoors |
| No. 815020 | + item No. |
| | |

| 156–196 cm | Recommended user height |
|---|---|
| 110 kg | Maximum weight of user |
| 5 k <u>ç</u> | Maximum load storage net |
| 81–100 cn | Handle height |
| Ergonomic, serveral holding position: (TOPRO ErgoBar | Handlebar |
| 7.5 cm | Grip distance driving brake |
| 37.5 cn | Width between center points of handles |
| 5 k <u>ç</u> | Maximum load tray |
| 68.3 cn | Height of tray position |
| 35 × 29 cm | Tray size |
| 78 cn | Turning circle |
| 100 cm | Maximum height |
| 55.5 cm | Maximum width |
| 63 cn | Maximum length |
| 11.2 cn | Diameter wheels |
| 2.4 cn | Width wheels |
| 25 cn | Length when folded |
| 55.5 cn | Width when folded |
| 88.5 cn | Height when folded |
| 25 × 55.5 × 88.5 cn | Dimensions folded $(I \times w \times h)$ |
| 6.9 kg | Product weight |
| Indoor | Area of use |
| 815500-119010 | Article No. |

Available in colour:



topromobility.com 19 18 topromobility.com

Technical data for Troja 5G

| Size | M (medium) | S (small) |
|--|--|--------------------|
| Recommended user height | 150–200 cm | 135–170 cm |
| Maximum weight of user | 150 kg | 125 kg |
| Maximum capacity shopping bag | 10 kg / 15 l | 10 kg / 15 l |
| Handle height | 78–100 cm | 69–86 cm |
| Handles | Ergonomic, serveral holding positions (TOPRO ErgoGrip) | |
| Grip distance driving brake | 7 cm | 7 cm |
| Width between center points of handles | 47 cm | 47 cm |
| Seating width | 45.5 cm | 45.5 cm |
| Seat height | 63 cm | 54 cm |
| Turning circle | 84 cm | 84 cm |
| Maximum height | 103.5 cm | 89.5 cm |
| Maximum width | 60 cm | 60 cm |
| Maximum length | 65 cm | 65 cm |
| Diameter wheels | 20 cm | 20 cm |
| Width wheels | 3.6 cm | 3.6 cm |
| Length when folded | 65 cm | 65 cm |
| Width when folded | 24 cm | 24 cm |
| Height when folded | 79.5 cm | 73 cm |
| Transport dimensions ($I \times w \times h$) (excl. wheels and shopping bag) | 65 × 24 × 73 cm | 65 × 24 × 66 cm |
| Transport weight (excl. wheels and shopping bag) | 5.6 kg | 5.4 kg |
| Product weight (incl. wheels and shopping bag) | 7.8 kg | 7.6 kg |
| Area of use | | Indoors & outdoors |
| Article No. | 814600 + item No. | 814601 + item No. |



Model M & S are available in on-trend matt coated colours:



Silver















Black

Technical data for Troja Original



| Size | M (medium) | S (small) | X (extra small) |
|---|-----------------|-----------------|-------------------|
| Recommended user height | 150–200 cm | 130–170 cm | 120–160 cm |
| Maximum weight of user | 165 kg | 165 kg | 125 kg |
| Maximum capacity basket | 10 kg / 20 l | 10 kg / 20 l | 5 kg / 11 l |
| Handle height | 78–100 cm | 69–86 cm | 64-79 cm |
| Handles | | | Anatomical |
| Grip distance driving brake | 6.5 cm | 6.5 cm | 6.5 cm |
| Vidth between center points of handles | 47 cm | 47 cm | 37 cm |
| Seating width | 45.5 cm | 45.5 cm | 36 cm |
| Seat height | 63 cm | 54 cm | 51 cm |
| Turning circle | 84 cm | 84 cm | 79 cm |
| Maximum height | 101 cm | 87 cm | 80 cm |
| Maximum width | 60 cm | 60 cm | 51 cm |
| Maximum length | 65 cm | 65 cm | 65 cm |
| Diameter wheels | 20 cm | 20 cm | 20 cm |
| Width wheels | 3.6 cm | 3.6 cm | 3.6 cm |
| Length when folded | 65 cm | 65 cm | 65 cm |
| Width when folded | 24 cm | 24 cm | 24 cm |
| Height when folded | 80 cm | 74 cm | 66 cm |
| Dimensions folded $(I \times w \times h)$ | 65 × 24 × 80 cm | 65 × 24 × 74 cm | 65 × 24 × 65 cm |
| Product weight (incl. wheels and basket) | 6.8 kg | 6.6 kg | 6.1 kg |
| Area of use | | | ndoors & outdoors |
| Article No. | 815400-110000 | 815410-111000 | 815420-112000 |

Model M, S and X are available in on-trend matt coated colour:



Silver

20 topromobility.com 21

Technical data for Olympos ATR



| Size | M (medium) | S (small) |
|--|--|-------------------|
| Recommended user height | 150–195 cm | 140–170 cm |
| Maximum weight of user | 150 kg | 150 kg |
| Maximum capacity shopping bag | 15 kg / 12 l | 15 kg / 10 l |
| Handle height | 76–95.5 cm | 69-86.5 cm |
| Handles | Ergonomic, serveral holding positions (TOPRO ErgoGrip) | |
| Grip distance driving brake | 7 cm | 7 cm |
| Width between center points of handles | 49 cm | 50 cm |
| Seating width | 47 cm | 47 cm |
| Seat height | 60 cm | 54 cm |
| Turning circle | 97 cm | 97 cm |
| Maximum height | 99 cm | 90 cm |
| Maximum length | 77 cm | 77 cm |
| Maximum width | 68 cm | 68 cm |
| Diameter wheels | 28 cm / 23 cm | 28 cm / 23 cm |
| Width wheels | 4.2 cm | 4.2 cm |
| Length when folded | 77 cm | 77 cm |
| Width when folded | 27 cm | 27 cm |
| Height when folded | 87.5 cm | 84 cm |
| Dimensions folded $(I \times w \times h)$ | 77 × 27 × 87.5 cm | 77 × 27 × 85 cm |
| Product weight (incl. wheels and shopping bag) | 8.9 kg | 8.7 kg |
| Area of use | | Outdoors |
| Article No. | 814300 + item No. | 814305 + item No. |
| | | |

available in on-trend matt coated colours:

Silver
Sand
Yellow
Wine red
Purple
Blue
Turquoise
Black

Model M & S are

The future is now
Choose a TOPRO rollator and help ensure a sustainable future.
Read more about our contribution on topromobility.com

Global responsibility

Sustainable materials and manufacturing

Clean energy

Sustainability in focus

We embrace the future by reducing our environmental footprint and contributing more to the circular economy, with implementing new accessory packaging made from materials that are easy to recycle.



22 topromobility.com 23

Premium rollators from Norway

Manufacturer:



TOPRO Industri AS Rambekkvegen 5 2816 Gjøvik Norway Tel. (+47) 61 13 46 00 info@topromobility.com www.topromobility.com





Designed and manufactured in Norway.

TOPRO vouches for top quality and offers a 7-year warranty.

Setup and variations in fixed accessories for the rollators may differ between countries.

Quality Assurance Standards

EN ISO 9001 – quality EN ISO 14001 – environment EN ISO 11199-2 – rollator EN ISO 10993-5 – biocompatibility

made in Norway



Dealer