

For an Exceptional Bathing Experience

The Bellavita NOVA offers superb comfort for an enhanced bathing experience, in a product which is very easy to install, very reliable and provides a safe and convenient way to support the continued use of your bath. Designed, developed and manufactured by our bathlift experts in Germany, the Bellavita NOVA offers several new exciting patented features.



optimum stability when rising and lowering across the full height range.



## **Key Features and Benefits;**

 Features a patented removable backrest which can be taken out by a carer when the bathlift is in its lowest position to maximise bath space, user comfort and water immersion.

6.5cm

A high seat

height of 48cm

enabling easy and

safe transfers in

and out of the

bath

- A patented dual rack ensures optimum stability when rising and lowering across the full height range, and excellent stability for transfers on and off the bathlift
- At the touch of a button, the Bellavita NOVA gently lowers users down to the bottom of the bath. After bathing and relaxation, the bathlift smoothly raises the user back to the top to enable a safe exit
- The backrest reclines to an incredible maximum angle of 50 degrees, operated easily by the hand control and enabling the user to relax in comfort
- A front cut out in the seat improves access and comfort for washing and personal hygiene

- Supports a maximum user weight of 170kg (26.7st)
- Very low seat height of 6.5cm (2.5") from bottom of bath for maximum water immersion, ideal for shallow baths
- A high seat height of 48cm (19")
   enabling easy and safe transfers in
   and out of the bath. No need for height
   adapters.
- An ergonomically designed backrest for lumbar support and added comfort
- At 9.5kg without the handset (9.9kg with the handset) the Bellavita NOVA is very lightweight making it easy to install and remove from the bath
- The Bellavita NOVA separates easily and quickly into two lightweight sections for ease of fitting/removal, storage and transportation
- The backrest includes an ergonomically designed handle to support removal of the backrest when in use and make installation and transportation easier and safer.



Lightweight, waterproof floating handset

- Contoured, slim-line profile gives a modern, aesthetically pleasing appearance to suit most bathrooms
- Four patented 3D suction feet to enhance fitting to the bath tub, maximising adhesion, stability and safety
- The front two are self-venting and flexible, enabling easier installation and removal of the bathlift
- Lightweight, waterproof and floatable handset with colour coded buttons is easily operated (Both the handset and motor are IP X7 rated)
- Hand control recharges quickly for convenience and multiple use environments
- An indicator light notifies users when the handset needs charging. The bathlift has a built in safety feature which prevents the bathlift from lowering down if there is not sufficient charge to return the bathlift seat back to the top of the bath

- Emergency stop button is incorporated into hand control for safety
- Lightweight plastic frame is easy to clean, with smooth surfaces and fits into most baths easily
- Can be used with or without covers; washable hygienic covers are easily fitted/removed and provide comfort when bathing (available in white or blue)
- Wide range of optional accessories; including covers, wider side flaps, turning and transfer aids, side flap deflectors, and a lap belt and chest strap.
- Shares a number of parts and accessories with the original Bellavita bathlift

Separates into two lightweight sections

Maximise bath space, user comfort and water immersion for a relaxing bathing experience

Patented
removable
backrest which can
be removed by
a carer





## Codes:

Code	Description	
470400100	Bellavita NOVA without covers	
470400252	Bellavita NOVA with classic white cover	
470400312	Bellavita NOVA with classic blue cover	








**Cover Set Classic Blue** 

## **Specifications:**

Overall Length (Reclined Backrest): <b>102cm</b>	Overall Width: <b>70cm</b>	
Seat Width: 38cm	Seat Depth: 52.5cm	
Backrest Height: 60cm	Backrest Width: 30-32cm	
Angle of Backrest Recline: 50°	Min. Seat Height: 6.5cm	
Max. Seat Height: 48cm	Width of Baseplate: 30cm	
Depth of Baseplate: 60cm	Product Weight (without handset): 9.5kg	
Product Weight (with handset): 9.9kg	Max. User Weight: 170kg	

## If you would like to view this product, please contact your local stockist.

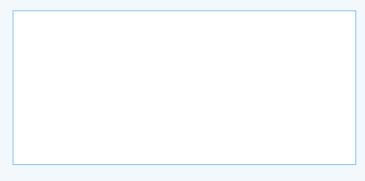



Code	Description
460800312	Cover Set Classic Blue
460800252	Cover Set Classic White
412100100	Vitaturn K Turning Aid
412100200	Vitaturn L Turning and Transfer Aid
460900227	5cm Side Flap Extension
460900228	7.5cm Side Flap Extension
460900229	10cm Side Flap Extension
400570	Side Flap Deflectors
460900220	Suction Cup set 4 x Grey
460900225	Suction Cup set 4 x Clear
460900150	4 Point Chest Strap
460900100	2 Point Chest Strap
410900100	Active Foam Cleaner





www.drivedevilbiss.co.uk



