

E-MED
Patient Trolley
Range



drive

 DeVilbiss[®]
HEALTHCARE



ENTRANCE



Introduction

Drive DeVilbiss Healthcare are one of the leading manufacturers and distributors of medical devices throughout the world.

The E-MED Range brings together a versatile and practical selection of patient trolleys, developed for applications in the acute care environment including A&E, hospital wards, day case surgery, x-ray and general patient transportation.

From the high performance E-MED Seatline for day surgery and critical care to the variable height E-MED 1400 & 1200 and the fixed height E-MED 1602 for general transfers. These high specification patient

trolleys have been designed and engineered to provide safe and comfortable surfaces, enabling diagnostic and treatment procedures in the acute care environment, reducing the need for surface to surface transfers.

Designed around flexibility, safety and performance, the E-MED patient trolleys offer a comprehensive range and choice, combining performance, safety and reliability.

Contents


E-MED Patient Trolley Range

- 06 E-MED Seatline
- 08 E-MED One Day EST
- 10 E-MED 1200 & 1400
- 12 E-MED 1500 Range
- 14 E-MED 1600
- 16 E-MED Accessories



- + +
E-MED
Patient Trolley
Range

E-MED Seatline /
E-MED One Day EST /
E-MED 1200 & 1400 /
E-MED 1500 Range /
E-MED 1600

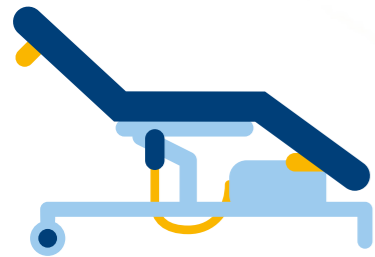

E-MED Seatline

“ Patient platform operated pneumatically or electrically ”

Power Drive option available for single finger manoeuvrability



Features & Benefits

The E-MED Seatline is the flagship model in our trolley range. The high specification Seatline is a multifunctional three-section trolley, suitable for a wide range of applications but developed specifically for day surgery and high dependency use.

Designed around optimum flexibility, ease of use, occupant comfort and safety, the Seatline is ideally suited to even the most critical applications. The structure of the trolley, which is endowed with a single central column for the raising adjustment, allows many x-ray devices (c-arm, for example) to be easily positioned under the backrest and leg rest sections. The result is brilliant: without moving the patient from the trolley, it is possible to operate x-ray exams of high quality and precision.

- ▶ Designed to allow insertion of a standard or digital x-ray cassette directly below the platform, facilitating high quality x-rays without moving the patient
- ▶ Variable height platform - descends to 65cm, rising to 99cm to facilitate treatment procedures in a day surgery environment
- ▶ Incorporates a 70° backrest incline and 40° leg rest incline, enabling positioning and angling for individual comfort and safety
- ▶ Adjustable safety side rails can be dropped to below the level of the platform using a simple motion assisted hand control
- ▶ 150mm antistatic braked castors ensure the Seatline is both manoeuvrable and smooth in use
- ▶ Mattress is breathable, water resistant, antistatic and antibacterial
- ▶ Central breaking system from any corner of the trolley
- ▶ Full range of accessories available

Comfortable true chair position trolley



Specifications

- Length: 215cm
- Width: 80cm
- Height (min): 65cm
- Height (max): 99cm
- Mattress Platform Length: 200cm
- Mattress Platform Width: 60cm

E-MED 2500
Product Code: EMED/2500
 (3 section, electric with power drive)

E-MED 2450
Product Code: EMED/2450
 (3 section, electric with manual drive)

E-MED 2400
Product Code: EMED/2400
 (3 section, pneumatic with manual drive)

 Lockable Handset (2500/2450)	 X-Ray Translucent	 Safe Working Load 180kg 28st
 3-Section	 2 Year Warranty	 Max User Weight 155kg 24st

E-MED One Day Extended Stay Trolley



“ Patient platform electrically operated ”

Features & Benefits


The E-MED One Day Extended Stay Trolley is an electrically operated four-section trolley suitable for a wide range of applications but developed specifically for extended stay use. This trolley has a minimum platform height of 40cm, the lowest of our range, and a maximum height of 83cm.

- ▶ Platform length extends a further 18cm to accommodate taller patients
- ▶ Flexible controls for the patients comfort and ease
- ▶ Mattress included is water resistant, antibacterial and multi-stretch. It also has a removable cover
- ▶ Central break from either side of the trolley
- ▶ 200mm antistatic castors as standard
- ▶ Full range of accessories available

Specifications

- Length: 199 - 217cm
- Width: 89.5cm
- Height (min): 40cm
- Height (max): 83cm
- Mattress Platform Length: 184 - 202cm
- Mattress Platform Width: 71cm

E-MED EST
Product Code: EMED/ONEDAY/EST
(4 section, electric)

 Lockable Handset	 260kg 40st	 4-Section
 2 Year Warranty	 225kg 35st	



Length extendable patient platform



E-MED 1200 & 1400



“
Patient platform pneumatically operated
”



Open architecture for easy c-arm access

Features & Benefits

The E-MED 1200 two-section and the E-MED 1400 four-section patient trolleys have been developed for applications in the acute care environment including A&E, hospital wards, day case surgery, x-ray and general patient transportation.

The trolleys operate with a smooth, adjustable height range of 56-88cm, easily adjusted using dual sided foot controls which can be placed at the head or foot of the trolley.



The special patient platform is totally radio transparent, and the open construction of the trolley beneath the platform enables x-ray and imaging equipment to be operated either above or below the full length of the patient.

- ▶ Minimal transfer gap with no lateral protrusions
- ▶ Full length x-ray cassette options
- ▶ 200mm antistatic castors as standard, 5th wheel option for easy movement through 360°
- ▶ Height adjustment pedals on both sides of trolley
- ▶ Patient documents/personal items spaces, oxygen tank holder as standard
- ▶ High quality 8cm anti-static, pressure reducing mattress
- ▶ Central braking system from any corner of the trolley
- ▶ Full range of accessories available



Specifications

- Length: 213cm
- Width: 83cm
- Height (min): 56cm
- Height (max): 88cm
- Mattress Platform Length: 194cm
- Mattress Platform Width: 62cm


E-MED 1200
Product Code: EMED/1200
(2 section, pneumatic)

E-MED 1400
Product Code: EMED/1400
(4 section, pneumatic)

E-MED 1500 Range



“ Patient platform pneumatically operated ”

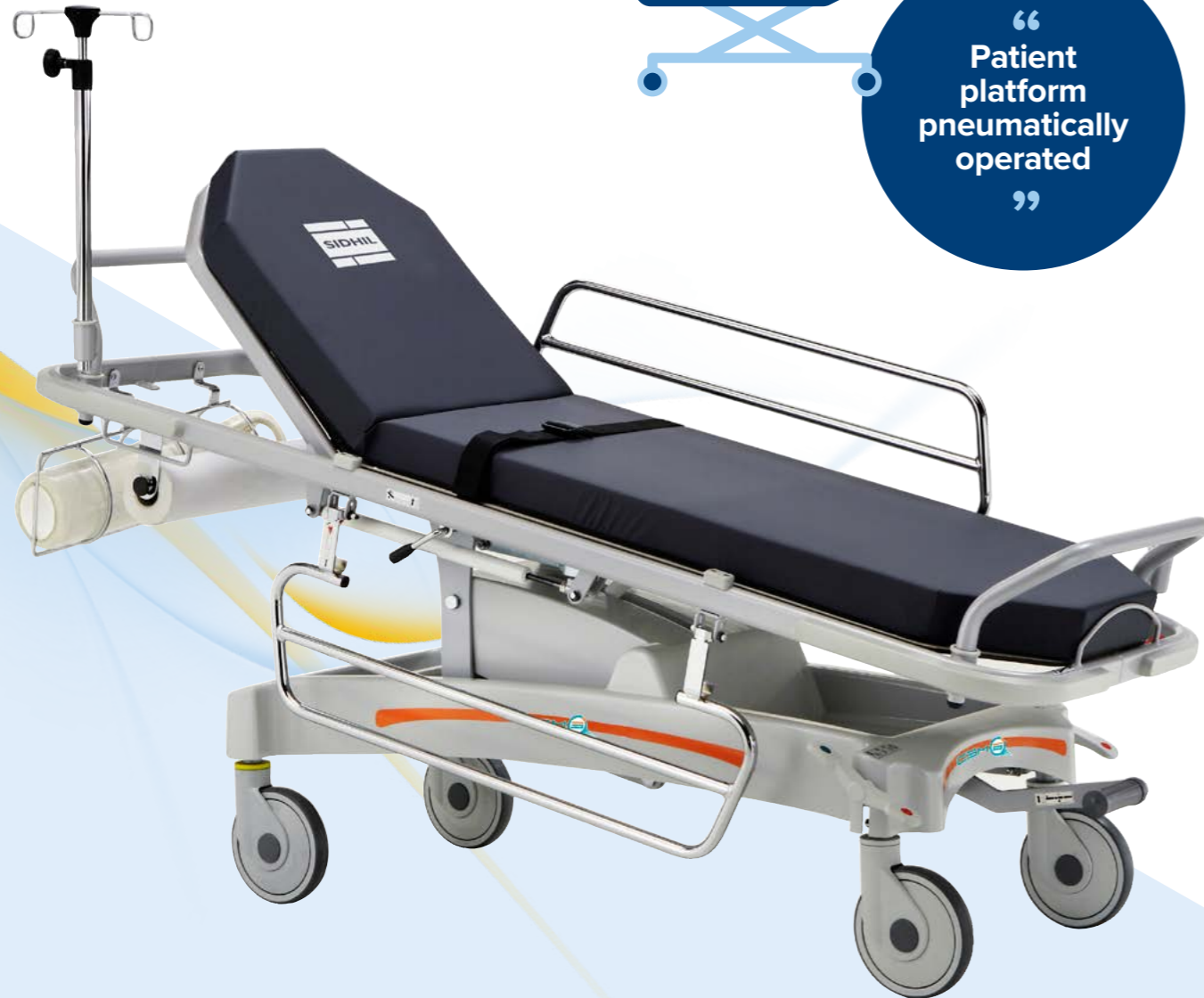

Features & Benefits

The E-MED 1500 Range are two and four section trolleys available with or without trendelenburg positioning. These trolleys have been developed for applications in the acute care environment including A&E, hospital wards, day case surgery, x-ray and general patient transportation.

The trolley is designed with a robust central lifting mechanism for height adjustment allowing x-ray devices (c-arm, for example) to be easily positioned under the backrest and leg rest sections.

- ▶ Two section trolley (backrest, pelvic/femoral/leg rest section)
- ▶ Four section trolley (backrest, pelvic/femoral/leg rest section with knee break)
- ▶ 200mm standard castors with optional antistatic addition, 5th wheel option for easy movement through 360°
- ▶ Radio transparent surface available as a factory fitted option
- ▶ A range of optional push handles and side rails available to ensure trolleys meet clinical needs
- ▶ For ease of use, adjustment levers are located on both sides of the trolley
- ▶ Full range of accessories available

Specifications

Length: 204cm

Width: 70cm

Height (min): 53cm

Height (max): 90cm

Mattress Platform Length: 193cm

Mattress Platform Width: 58cm

E-MED 1500

Product Code: EMED/1500
(2-section, pneumatic, no trendelenburg)

E-MED 1504

Product Code: EMED/1504
(2-section, pneumatic, with trendelenburg)

E-MED 1510

Product Code: EMED/1510
(4-section, pneumatic, no trendelenburg)

E-MED 1512

Product Code: EMED/1512
(4-section, pneumatic, with Trendelenburg)



X-Ray Translucent
(Factory Fitted Option)



Safe Working Load
300kg
47st



2-Section
(1500/1504)



4-Section
(1510/1512)



2 Year Warranty



Max User Weight
265kg
42st



Optional 4-section profiling platform



E-MED 1600



“ Patient platform pneumatically operated ”



Features & Benefits

The E-MED 1600 range are lightweight, fixed height trolleys featuring two separate platform sections. The backrest section has multiple positions and can be adjusted manually up to 80° for patient comfort.

- ▶ Fixed height 2 section patient trolley of 78cm
- ▶ Backrest with multiple positions for patient comfort can be positioned up to 80°
- ▶ Lightweight frame with side bumpers enables easy trolley movement
- ▶ 200mm castors as standard
- ▶ A range of optional push handles and push bars available
- ▶ Full range of accessories available



Specifications

Length: 204cm

Width: 70cm

Trolley Height: 78cm

Mattress Platform Length: 193cm

Mattress Platform Width: 58cm

E-MED 1600

Product Code: EMED/1600
(2-section, pneumatic, independent brake)

E-MED 1602

Product Code: EMED/1602
(2-section, pneumatic, constant brake)



Safe Working Load



2-Section



2 Year Warranty



Max User Weight



Constant brake option available

Compass and folding side rails available



Contact Us



For further information, please contact:

Drive DeVilbiss Healthcare Ltd
Heathfield Lane, Birkenshaw
West Yorkshire, BD11 2HW
T: 0845 0600 333
E: info@drivedevilbiss.co.uk
W: www.drivedevilbiss.co.uk

Enhancing the
Quality of Life of the
People We Touch.



drive

 **DeVilbiss[®]**
HEALTHCARE

www.drivedevilbiss.co.uk

Leaflet Code: LL430

All sizes and weights stated are nominal. All errors and omissions excepted. Due to the limitations of photography and the printing process colours shown may not be 100% accurate. All equipment should be used in conjunction with clinical judgement and assessments dependent on individual patient morphology. Essential nursing care is pivotal in providing the optimum pressure ulcer prevention plan. *The risk factors stated in this brochure are used in conjunction with the NICE Guidelines